

HOOF CARE PRODUCTS AVAILABLE AT AN FPD DEALER NEAR YOU

# IT'S ALL ABOUT THE HEALTH OF THE FOOT




KERCKHAERT

[WWW.FARRIERPRODUCTS.COM](http://WWW.FARRIERPRODUCTS.COM)



## STEEL COMFORT SPORT

Kerckhaert front shoe designed for ease of breakover and good support. The beveled inside web of shoe promotes cleanout. Improves coordination of movement and optimized balance with less pressure on the tendons. It offers a 5-nail pattern for extra nailing options and a V-Crease for more secure nail fit to hold shoes better. It is already marked to drill for traction devices in the heel area.

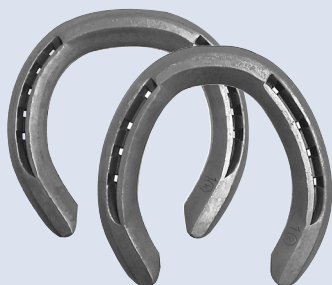


## CLASSIC ROLLER

While a half round style, the Kerckhaert Classic Roller has more ground surface than a true half round. The Classic Roller is available in quarter clipped or toe clipped front, and side clipped hind.

## ALUMINUM COMFORT

The Kerckhaert Aluminum Comfort has a rolled edge built into the toe and sides to enhance breakover and allow for easier lateral movement. Clips provide additional security.

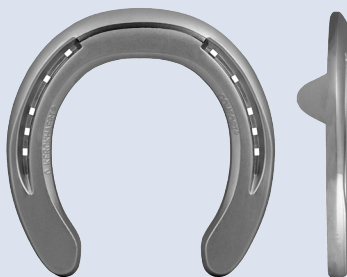

# PLEXUS PERFORMANCE PADS



Made in the U.S., Austin Edens' Plexus Pads feature a patent pending, common sense design. Plexus Pads are designed to fit both fronts and hinds. Shape and design allows for easy fit for front or hind for more hoof sizes, reducing need for more inventory. High quality injection molded material allows for resets.



## OPEN TOE CAUDAL SUPPORT BENEFITS

- Unique diamond-shaped holes provide greater locking strength for FootPro™ DIM and Vettec Equi-Pak
- Frog support is tapered at end of pad to reduce the need to trim or grind pads, or grind heels of shoes
- Open Toe design avoids trapping of foreign materials in front half of foot



**FLAT**



**DEGREE**




## FULL PAD BENEFITS

- The Full Plexus pads are a hybrid design, size S-M, to minimize the number of sizes you need to carry.
- The one size will easily fit shoe sizes 00 through 3



**FULL  
FLAT**

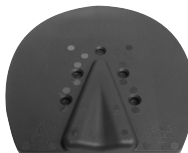


**FULL  
DEGREE**





## DIAMOND MAX GLIDE RASP

The innovative and unique manufacturing process of the Diamond Max Glide creates a high-quality rasp and file; 10% wider than the standard rasp. The rasp side features a 9-tooth pattern with slightly smaller teeth, and a file side with Chip Breaker.

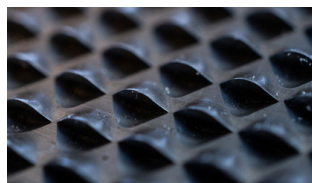

## PLEXUS CS GEL

Plexus CS Gel has been developed by Austin Edens to help treat white line and other bacterial/fungal issues in the horse's foot. The gel formula contains copper sulfate and vinegar in a ratio that makes it an effective treatment and preventative product.



# FOOTPRO™ HOOFCARE



## FOOTPRO™ CAST

FootPro Cast is a resin impregnated fiberglass mesh casting material that is an affordable and user-friendly option for the treatment of hoof problems where support and protection are needed. FootPro Cast provides protection from wear, helps strengthen and stabilize the hoof capsule and is excellent for therapeutic applications. It is wrapped dry with the help of adhesive holding the first layer, then wetted to activate the resin and cure it.



## FOOTPRO™ DIM

FootPro DIM Impression Material is a non-toxic, non-allergenic, dental grade silicone, hand mixed to a moldable consistency and applied to the bottom of the equine foot with a pad to provide either a soft cushion (DIM20), or a firm supportive material (DIM40) when needed. Each box contains a total of either 4 lbs. or 2 lbs. (Equal Parts Catalyst and Base).



## FOOTPRO™ CS+ HOOFCARE

The FootPro CS PLUS® Hoof Treatment is a clay based hoof treatment with copper sulfate, iodine, venice turpentine and tea tree oil. The product is used to treat hoof separations or surface pockets in the sole and frog where bacteria may be present. Packed in a convenient 10 oz. container, the material is easily spread and leaves no stain on your fingers.



# VETTEC HOOFCARE

**Vettec**  
ANIMAL HEALTH



## EQUI-PAK

Equi-Pak is a liquid urethane that sets in 30 seconds and is used as a pour-in pad to provide comfort, protection and therapeutic support. It bonds directly to the sole and frog and can be used by itself as an instant pad, or as a packing underneath a pad. Equi-Pak limits exfoliation, promotes sole depth and makes the bottom of the hoof impervious to injury from sharp objects.



## EQUI-PAK SOFT

Equi-Pak Soft is a liquid urethane that sets in 40 seconds and is dispensed under a pad to provide comfort, protection and therapeutic support. It bonds directly to the sole, frog and pad to keep out debris.



## EQUI-PAK CS

Equi-Pak CS is a copper sulfate infused liquid urethane that sets in 60 seconds and is used as a pour-in pad to provide comfort, protection and therapeutic support for equine feet. It bonds directly to the sole and frog and can be used by itself as an instant pad, or as a packing underneath a pad. Equi-Pak CS limits exfoliation, promotes sole depth and makes the bottom of the hoof impervious to injury from sharp objects.



## SIL-PAK

Sil-Pak is a non-bonding, skin-safe silicone packing material dispensed under a pad to seal, protect and cushion the equine foot. It totally fills and seals all void areas and prevents debris migration. Silpak applications are weight bearing in 2 minutes.



## SIL-PAK SOFT

Sil-Pak Soft has all the features of Sil-Pak but with a shore hardness of 15 (Sil-Pak is approximately 35). Sil-Pak Soft applications are weight bearing in 2 minutes.



## ADHERE

A fast-setting urethane adhesive designed to bond steel, aluminum and urethane shoes to the hoof and to fabricate large and small hoof wall repairs. It sets in 60 seconds, is weight bearing after 3.5 minutes and is fully cured in 6 minutes. Perfect for most every hoof repair application, Adhere sets hard with a strong bond, yet maintains the flexibility necessary at the quarters and heels to allow the natural movement of the hoof. Black color dresses to gray.



## SUPER FAST

Super Fast is a fast setting, highly durable urethane formula that is ideal for making custom shoes directly on the hoof, as well as medial, lateral and toe extensions on foals. It sets in 30 seconds and is weight bearing in 2.5 minutes. Custom shoes made from Super Fast, when properly applied, should last 5 to 6 weeks or longer. Super Fast can also be used to repair hoof wall defects, externally modify shoes and glue – on wood and urethane clogs.



## SOLE-GUARD

Sole-Guard is a black liquid urethane that can be used in a variety of applications to produce an extremely durable pad on shod and unshod equine feet. It sets in 30 seconds, bonds to the sole and frog and can be used as a thin protective pad, or a thicker, supportive pad.



# HOOF CARE ACCESSORIES

## VETTEC 210CC DUAL SPEED DISPENSER

Vettec now has a 210cc Dispenser with a 2-speed feature. One setting of the knob is the same as the normal dispenser and the second setting dispenses at a much faster rate, putting out more material with less pumps.



**FOOTPRO 210CC  
DISPENSING GUN**



**DIAMOND 210CC  
APPLICATOR GUN**



**210CC  
MIXING TIPS**



## FOOTPRO™ HOOF BUFFER ATTACHMENT


Attaches to cordless or electric drill.  
2" x 4" sleeves in 60, 80 and 100 Grit.

## DEPLANO WOODEN SHOES

Typically used in treatment of chronic and sometimes acute laminitis cases. This is a complementary product for the FootPro Cast, as most wooden shoes or clog applications are done with a cast material to secure the shoe.





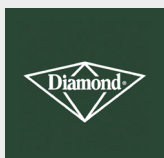

## DIAMOND FROG DRESSING

Diamond Frog Dressing is easy to apply and control thanks to its gel formulation. Treatment works against fungus and bacteria for maximum effectiveness against thrush causing organisms. Contains no carcinogens. The blue color shows where product has been applied, and when it needs to be reapplied. The flexible tip applicator gets treatment into cracks, nail holes, and defects.



## ALLEN'S BLUE POWDER

Allen's Blue Powder is a medicated thrush preventative, developed by farrier Todd Allen. The active ingredient is copper sulfate and the powder is sprinkled on the frog and sole area of the foot as a preventative measure. Also used as a treatment for thrush.

## HAWTHORNE SOLE PACK

Relieves dry, hard, sore hooves, while combating bacterial and fungal infections. It also treats white line disease. Extremely effective in maintaining the natural pliability of the hoof. Fast-acting, quick-penetrating and formulated with soothing natural ingredients. It's unique ability to bring the hoof back to a natural state makes it effective whether you're faced with a hoof that is too soft or a hoof that is too hard. Using Hawthorne Sole Pack a few days prior to trimming will bring the hoof to a pliable state that will ease the trimming process.





# HOOF STANDS & ACCESSORIES

## BASES, CRADLES & POSTS

Hoofjack bases are available in standard, medium and mini with a full line of complementary posts, cradles, dental attachments and replacement parts.



## HOOFJACK TABLE

12"x12" aluminum table fits in any Hoofjack base. Designed to keep items like Vettec and other hoof repair or glue-on products within reach when working under the horse.

## FOOTPRO™ LIMB WITH DIAMOND TRAINING HOOF

Fits on any Hoofjack base and includes the post. The Diamond Training Capsule (optional) mounts to the limb. The limb is made from Hickory Wood for durability and strength. This is an excellent teaching aid and can be mounted for hoof surface or work on outer wall.



## HOT FIT PLATE

Designed to fit one section of the Hoofjack base. This plate gives farriers who hot-fit a place to set the hot shoe while working on the hoof.



## MINI CRADLE WITH STANDARD HEIGHT POST

When rasping or nailing, the cradle is not in the way, but still provides excellent support for the hoof. Standard 11" post allows for taller adjustment.





# Professional Farriers Demand THE BEST

## FPD DOES TOO

*"Thank you for all you do to support us farriers. FPD, Inc. has gone above and beyond to help my business run smoothly and also inspire the young bloods that are interested in the industry."*  
- Mike Talbot

### WHAT'S NEW

#### PLEXUS PERFORMANCE PADS



#### KERCKHAERT TREK



#### KERCKHAERT PRIME FIT



#### PLEXUS CS GEL



#### FOOTPRO SOFT DIM20 IMPRESSION MATERIAL



#### FOOTPRO CARBIDE STUDS



#### LIBERTY SECURE HARD TOP NAILS



FIND AN FPD DEALER NEAR YOU CARRYING THESE EXCITING NEW PRODUCTS  
AT [FARRIERPRODUCTS.COM/LOCATIONS](http://FARRIERPRODUCTS.COM/LOCATIONS) OR CALL (800) 468-2879.



@farrierproductdistribution



@fpdinc



@farrierproducts





FARRIER PRODUCT DISTRIBUTION, INC. IS PROUD TO PROVIDE THE FINEST PRODUCTS AVAILABLE IN THE FARRIER INDUSTRY. INFORMATION ON HORSESHOES, NAILS, TOOLS AND RELATED EQUINE PRODUCTS CAN BE FOUND AT [FARRIERPRODUCTS.COM](http://FARRIERPRODUCTS.COM). IN ADDITION TO QUALITY PRODUCTS, WE ARE HAPPY TO PROVIDE YOU WITH PRACTICAL IDEAS AND SOLUTIONS FOR EVERYDAY APPLICATION – AVAILABLE UNDER THE NATURAL ANGLE AND FIELD GUIDE SECTIONS OF THE SITE.

PROFESSIONAL FARRIERS ACROSS NORTH AMERICA HAVE DISCOVERED THE VALUE IN USING PRODUCTS AVAILABLE THROUGH FPD DEALERS.

**FIND AN FPD DEALER NEAR YOU BY VISITING  
[FARRIERPRODUCTS.COM/LOCATIONS](http://FARRIERPRODUCTS.COM/LOCATIONS).**



@fpdinc



@farrierproductdistribution



@farrierproducts

